




RIKI GLOBAL
performance through innovation
Authorized Distributor

HPLC Columns

Paclitaxel

ACE C18-HL 5 μ m 250 x 4.6 mm

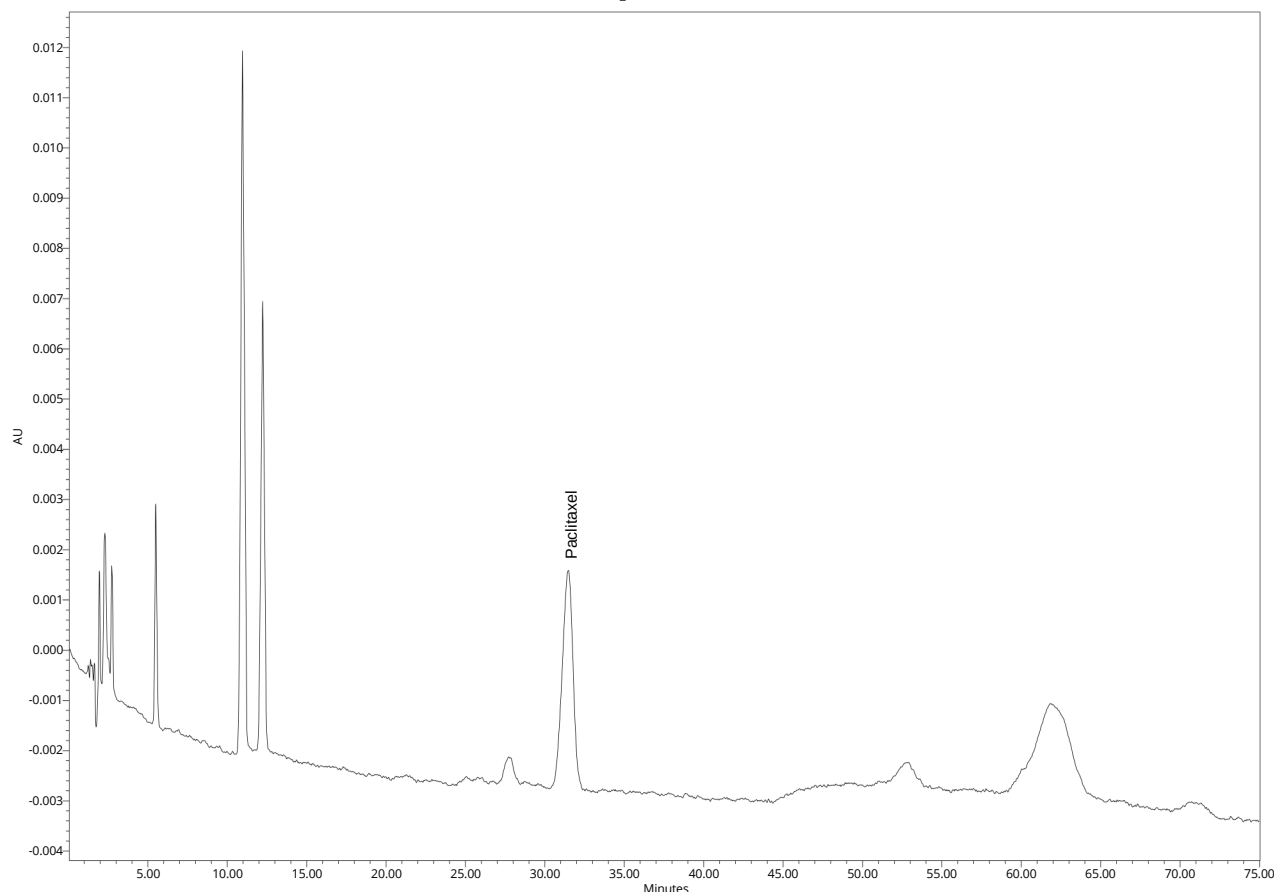

Test Condition

Column: ACE C18-HL 5 μ m 250 x 4.6 mm

Part No.: ACE-321-2546

Flow Rate: 1.2 mL/min

Injection volume: 10 μ L

Detection: PDA 227 nm

Mobile Phase: Water: Acetonitrile (95:5 V/V)

Column Temperature: 35°C

Sample Temperature: 25°C

For more information, Contact us at

✉ info@rikiglobal.com

🌐 www.rikiglobal.com

☎ +91 99082 46695

☎ +91-40-35172288/89